

ERA®-RV

Micro ERA®

Introducing the smallest extracoronal resilient attachment in the world!


Micro ERA®



The **ERA®-Reduced Vertical male (ERA®-RV)** has 0.4mm of vertical resiliency and universal joint hinging. The **Micro ERA® male** has the same resiliency and hinging, but needs 0.5mm less vertical space and has a diameter of almost 1.0mm less, making it the smallest extracoronal resilient attachment in the world. And yet, testing shows no appreciable loss of retention. The **ERA®-RV** and **Micro ERA®** males also have a projection which contacts the abutment crown above the female eyelet. This resists vertical displacement of the partial denture's distal extension saddle.

- Smallest extracoronal semi-precision attachment in the world that incorporates universal hinging with vertical resiliency
- Female is waxed into the pattern of the abutment restoration for a one-piece casting
- Black fabrication male with built-in spacer
- Six color coded males for six levels of retention (lightest to strongest): white, orange, blue, grey, yellow, red
- Nylon male greatly reduces the rate of female wear
- Optional Micro ERA® Metal Jacket holds the attachment male in the denture base and is sold pre-loaded with a black processing male
- Males easily changed without use of autopolymerizing acrylic

Female retention of the original white male, indicative of both male and female wear.

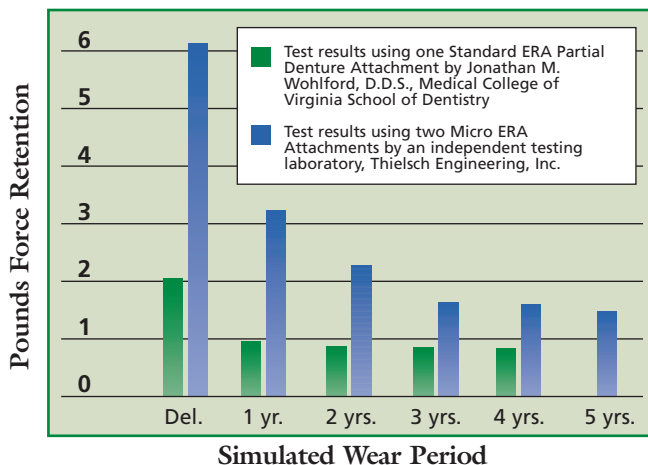

Micro ERA® Starter Kit

#811011

Kit contains: 2 complete attachments, 2 processing jigs, 1 core cutter bur, 1 seating tool, and 1 paralleling mandrel



Sterngold™

Order online at www.sterngold.com

Customer Service: 800-243-9942 • Toll Free Fax: 800-531-2685
 23 Frank Mossberg Drive, P.O. Box 2967 • Attleboro, MA 02703-0967

New Product

Special offer available for a limited time only.