

AHP-85MB & AHP-85MB-C HANDPIECES

IRRIGATION TUBING SETUP GUIDE

This guide provides step-by-step instructions for installing internal and external irrigation tubing onto the AHP-85MB and AHP-85MB-C Mont Blanc® contra angle handpieces.

(NOTE: AHP-85MB handpiece shown.)

Step 1: Remove Mont Blanc® handpiece (see Figure 1, Item #1) from box and install onto dental motor (consult motor operator's instructions). Insert rotating instrument into head of handpiece until "click" is heard.




Step 2: Preassemble Items #2, #3, #4, and #5 into a subassembly, as shown in Figure 2.

Step 3: Mount subassembly from Step 2 onto handpiece/motor, as shown in Figure 3.

Step 4: Connect tubing components with tubing from irrigation pump, as shown in Figure 3.

Single-Irrigation Setup: To set up irrigation for either internal only, or external only, connect the tubing w/small ID directly to the desired internal or external port — do not use the 'Y' connector. See Figure 4.

NOTE: When tubing is set up for simultaneous internal/external irrigation (Figure 3), the irrigation fluid will favor the line of least resistance, so the bulk of the fluid will flow to the external port.

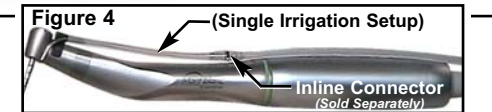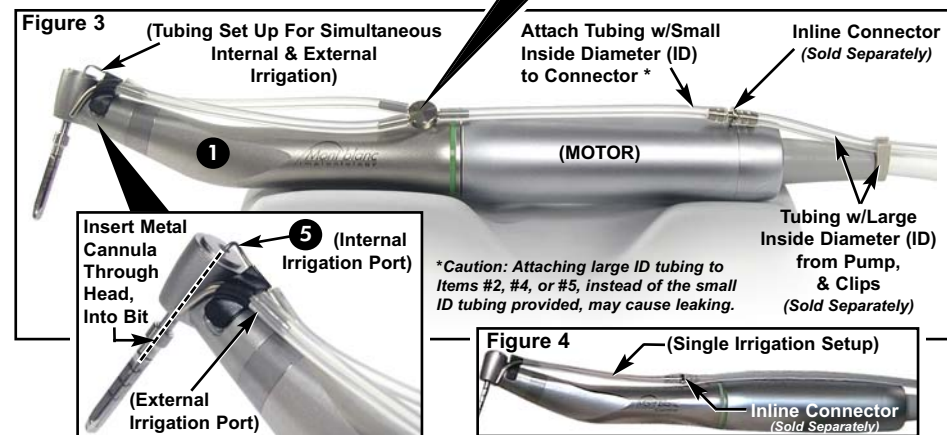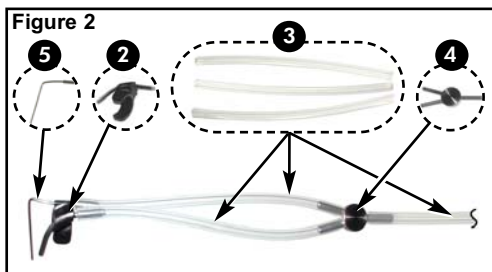

AHP-85MB & AHP-85MB-C HANDPIECES

IRRIGATION TUBING SETUP GUIDE

This guide provides step-by-step instructions for installing internal and external irrigation tubing onto the AHP-85MB and AHP-85MB-C Mont Blanc® contra angle handpieces.

(NOTE: AHP-85MB handpiece shown.)

Step 1: Remove Mont Blanc® handpiece (see Figure 1, Item #1) from box and install onto dental motor (consult motor operator's instructions). Insert rotating instrument into head of handpiece until "click" is heard.




Step 2: Preassemble Items #2, #3, #4, and #5 into a subassembly, as shown in Figure 2.

Step 3: Mount subassembly from Step 2 onto handpiece/motor, as shown in Figure 3.

Step 4: Connect tubing components with tubing from irrigation pump, as shown in Figure 3.

Single-Irrigation Setup: To set up irrigation for either internal only, or external only, connect the tubing w/small ID directly to the desired internal or external port — do not use the 'Y' connector. See Figure 4.

NOTE: When tubing is set up for simultaneous internal/external irrigation (Figure 3), the irrigation fluid will favor the line of least resistance, so the bulk of the fluid will flow to the external port.

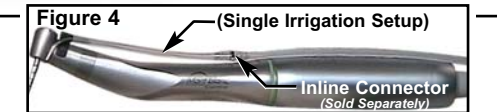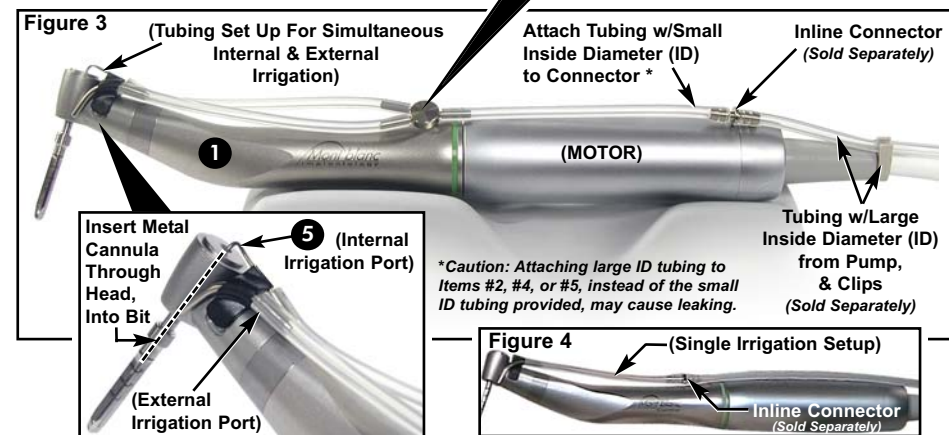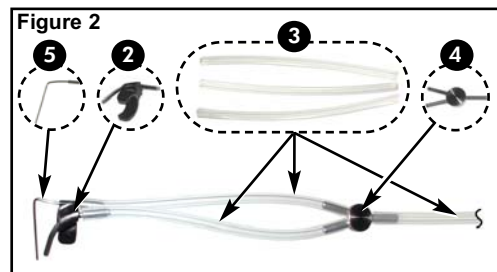

Sterngold™

Sterngold Dental
23 Frank Mossberg Drive • Attleboro, MA 02703
508-226-5660 • 800-243-9942
Fax: 508-226-5473 • www.sterngold.com



Sterngold™

Sterngold Dental
23 Frank Mossberg Drive • Attleboro, MA 02703
508-226-5660 • 800-243-9942
Fax: 508-226-5473 • www.sterngold.com