
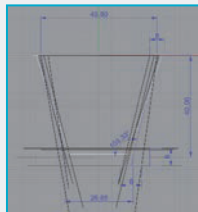


THE NEW INDIVIDUAL STANDARD

BD Load®

- The Baltic Denture System is based on the development of standardized set-up techniques, which are digitally stored in a software.
- In physical form, these have been implemented in a blank - the BD Load®.
- The teeth are already polymerized into the BD Load®. The design and milling process only serves to adapt the denture base.
- Milling out tooth compartments and bonding in denture teeth are not required, hence making the process economical and efficient.

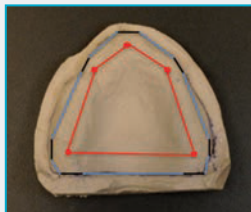

CHOOSE YOUR INDIVIDUAL STANDARD

BD Load® ECO

- 1 layer anteriors
- 1 layer posteriors

BD Load® Standard

- 4 layer anteriors
- 3 layer posteriors

- Sizes: S, M, L
- Jaw widths: narrow, medium and wide
- Set-up: either to the 1st or 2nd molar
- Shades: 16 V-Classic + 4 Bleach shades



OUR MASTERPIECE

- Lingualized occlusion according to Prof. Dr. Gerber with unilateral support
- Autonomous masticatory stability through tooth-to-tooth relationship
- Modern functional design and aesthetics
- Gentler centric through the mortar-pestle principle



Distributed by:

**DIGITAL
SOLUTIONS**



23 Frank Mossberg Drive
Attleboro, MA 02703
Tel: 1-800-243-9942
Fax: 1-508-226-5473
www.sterngold.com